


Offset Acetabular Reamer Handle

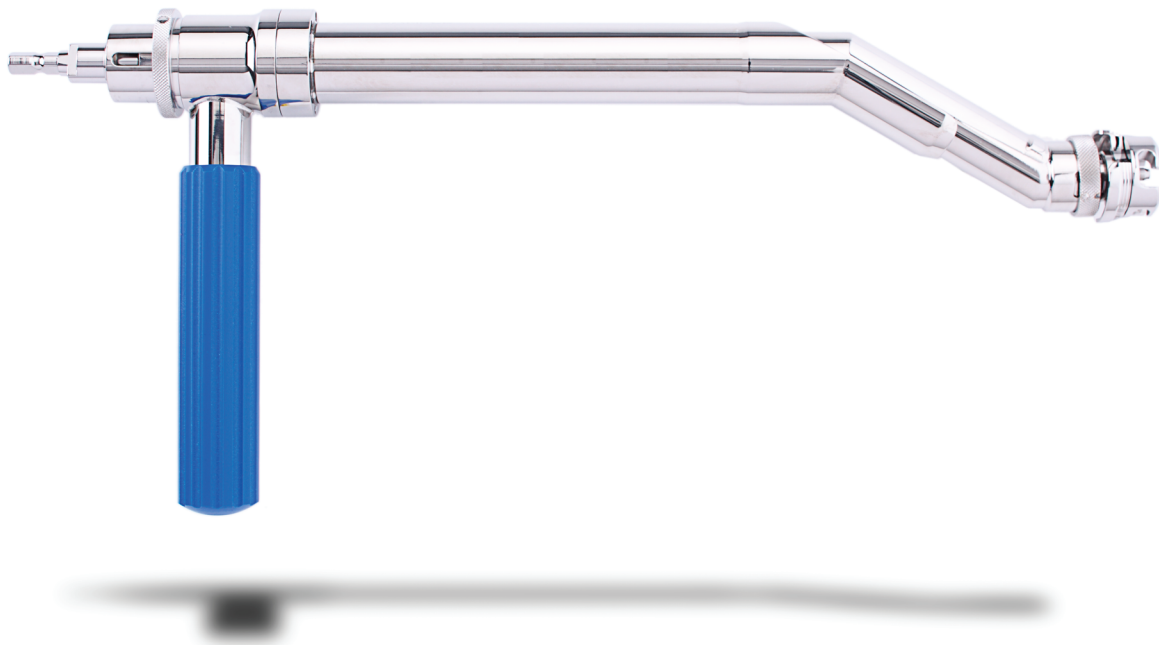

Inspired to facilitate less traumatic surgical procedures, the Offset Acetabular Reamer Handle incorporates design features that enable correct orientation reaming, while minimizing soft tissue damage.

RELIABLE PERFORMANCE

- The cross design features a simple and stable connection of the reamer with no movement, providing quick assembly and disassembly.
- The handle body separates into two halves to facilitate drive chain removal for optimum cleaning.

CUSTOMIZED PERFORMANCE

- A 30° offset of the reamer connection from the main handle allows the reamer to be correctly positioned in the acetabulum with the femur still occupying the space.
- The side handle adjusts for different angular orientations, providing multiple customized surgical approaches.

Patent Pending

Greatbatch Medical Headquarters
10000 Wehrle Drive
Clarence, NY 14031
U.S.A.
tel +1 716-759-5600
fax +1 716-759-2562

US Customer Service
4532 Park 30 Drive
Columbia City, IN 46725
U.S.A.
tel +1 260-244-6300
fax +1 260-244-6345

Europe & Asia Customer Service
L'Echelette 7
2534 Orvin
Switzerland
tel +41 32 358 01 11
fax +41 32 358 01 00

orthopaedicsna@greatbatchmedical.com

orthopaedicseurope@greatbatchmedical.com

www.greatbatchmedical.com